



MED SUPPLY LIFE

Handheld Ultrasound
Authorized Distributor

DOUBLE PROBE ULTRASOUND

AFFORDABLE
VISIBLE - FAST - EASY

MED SUPPLY LIFE CO., LTD.


Sales | Import | Installation | Maintenance of Medical Imaging Equipment

Model: D2CP - Handheld Double Probe Ultrasound

Tel / WhatsApp: +66 84-199-9893

LINE Official: @777dzaoj

E-mail: msupplylife@gmail.com

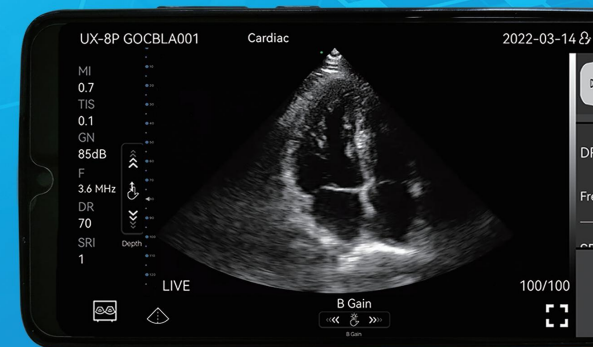

YOU CAN USE IT ANYTIME & ANYWHERE



APPLICATION

Application Sites: Cardiac, Superficial, Abdomen, etc.

Application Directions:

- Medical students training in college: Visual understanding of human body structure or ultrasound image study.
- Private clinic: Cardiac, Superficial, Abdomen, etc check.
- Preliminary examination: Home care, family planning.



FEATURES

- Highest cost-effective: Students can afford it, just like buying a cell phone.
- One scan whole body: D2CP composed of phased and linear array, it can scan from abdomen to superficial, especially heart.
- Small, light, built in battery: You can easily use it anywhere.

SPECIFICATION



- Scan Angle/Length: Phased 80, Linear 25mm
- Frequency: Phased probe 2.2/3.6MHz, Linear probe 7.5/10MHz
- Scan Mode: B, B/M, Color, PW
- Display Depth: Phased 90~190mm, Linear 20~100mm
- Power: Built-in lithium battery, wireless or usb charge 1.5h, work 1.5h
- Working System: iOS, Android, Windows
- WiFi Type: 802.11n/2.4G/5G
- Dimension: 132x55x20mm
- Weight: 150g



Innovative design

Linear/Convex/Phased Array 3 in 1, Linear/Phased Array, Pocket Ultrasound, Can be used in cardiac examination.



WIRELESS CHARGE

Easy and connect with a wireless power bank.



WiFi & USB CONNECT

Both USB and WIFI connected to mobile phones, tablets.



How it works?

- "Wireless ultrasound probe" is a mini ultrasound scanner that without screen. We made the main unit condensed into a small circuit board built-in the probe, and showing image in smart phone/tablet through Wifi transferring. - Image transferring through internal wifi from probe, no need external Wifi signal.

HIGH QUALITY IMAGES

See The Difference Great Image Quality Makes

